

TRIAL LENS SET & TRIAL FRAME

TLS-SERIES

TLS FULL DIAMETER



TLS DELUXE

TRIAL FRAME



D ELUXE TRIAL LENS SET WITH CORRECTED

Topcon knows what many practitioners look for in a trial lens set: The highest degree of accuracy to determine visual acuity, and convenience of use. The Topcon Deluxe Trial Lens Set offers a complete selection of lens powers, featuring the optical superiority of a special corrected curve design that compensates for thickness and air space between combined lenses. This allows the practitioner to read the exact correction directly from the frame.

The Topcon Deluxe Trial Lens Set also offers a variety of features for easy and efficient operation, including:

- Lenses marked only on one side to ensure correct positioning in the trial frame.
- Sturdy lightweight aluminum rings that are color-coded and marked for easy identification: red for minus, black for plus and blue-gray for prisms and accessories.
- All lenses are cemented as well as firmly clamped to ensure precision in the alignment.

The Corrective Curve Trial Lens Set comes in a sturdy holding tray with beveled slots for fast and easy insertion. The Deluxe Set includes minus and plus spherical lenses, prisms, other accessory lenses and various cylinder power lenses.



CARRYING CASE
optional

FULL DIAMETER TRIAL LENS SET

The Topcon Full Diameter Trial Lens Sets are a cost-effective combination of spheres, cylinders and accessory lenses. These economical bi-convex and bi-concave lenses are mounted in metal rims that have the lens power engraved in both sides of the handle. The full diameter lens provides a wide field of view during refraction. These lenses can be used with instruments that require a patient's correction such as a visual field. The Topcon Full Diameter Trial Lens Sets are supplied with a handsome velvet lined vinyl carrying case.

SERIES TLS SERIES

CURVE LENSES



Carrying Case is standard with the Full Diameter Trial Lens Set

Deluxe Trial Lens Set

- Spheres	Prisms	Accessories	+ Spheres
	20.00*	Red Maddox Rod*	20.00*
	12.00*		16.00* 18.00*
	8.00*		14.00* 15.00*
	5.00*	Pin Hole*	13.00* 13.50*
	4.00	Stenopaic Slit*	12.50
	3.00	Red Lens*	12.00
20.00*	2.00		11.50
16.00*	1.50	±0.50	11.00
13.00*	1.00	Cross	10.50
12.00	.50	Cylinders	10.00
11.00		Occluder	9.50
10.00		0.00 Sphere	9.00
9.00			8.50
8.00			8.00
7.50			7.50
7.00			7.00
6.50			6.50
6.00			6.00
5.50			5.50
5.00			5.00
4.50			4.50
4.00			4.00
3.75			3.75
3.50			3.50
3.25			3.25
3.00			3.00
2.75			2.75
2.50			2.50
2.25			2.25
2.00			2.00
1.75			1.75
1.50			1.50
1.25			1.25
1.00			1.00
.75			.75
.62			.62
.50			.50
.37			.37
.25			.25
.12			.12

*Individual Lenses Total Lenses : 271 pcs.

Full Diameter Trial Lens Set

- Spheres	Prisms	Auxiliary	+ Spheres
20.00	10.00*	.50[+/-]Jackson	20.00
18.00	8.00*	Cross Cylinder*	18.00
16.00	6.00*	Red-Free Lens*	16.00
15.00	5.00*	Green-Free Lens*	15.00
14.00	4.00	Cross Hair lens*	14.00
13.00	3.00	Plano Lens*	13.00
12.00	2.00	1.00mm Pinhole*	12.00
11.00	1.00	1.50mm Pinhole*	11.00
10.00	.50	1.00mm Slit*	10.00
9.00		Maddox Lens*	9.00
8.00		Occluder*	8.00
7.50			7.50
7.00			7.00
6.50			6.50
6.00			6.00
5.75			5.75
5.50			5.50
5.25			5.25
5.00			5.00
4.75			4.75
4.50			4.50
4.25			4.25
4.00			4.00
3.75			3.75
3.50			3.50
3.25			3.25
3.00			3.00
2.75			2.75
2.50			2.50
2.25			2.25
2.00			2.00
1.75			1.75
1.50			1.50
1.25			1.25
1.00			1.00
.75			.75
.50			.50
.25			.25
.12			.12

*Individual Lenses Total Lenses : 260 pcs.

T RIAL FRAME


The Topcon all-metal Trial Frame is constructed for durability and long life, while being lightweight and comfortable to wear.

The practical screw-thread design on all fittings allows for easy replacement of springs, templates and other components when the frame requires service.

The leaf spring loaded temples hold the frame snugly without creating any discomfort for the patient.

For a perfect fit and accurate settings, the TOPCON Trial Frame offers a vast array of features such as:

- Horizontal and vertical bridge adjustments.
- An adjustable saddle bridge assures a comfortable nose fit.
- Separate PD adjustment for each eye to compensate for asymmetrical factors in facial structure.
- Individual adjustments for temple length and angle.
- Adjustment for rotating cylinders to correct the axis.
- Scales with large, easy-to-read numerals.
- The spacing of the lens holders ensures accurate additive readings when the TOPCON Corrected Curve Additive Design lenses are used in the frame.

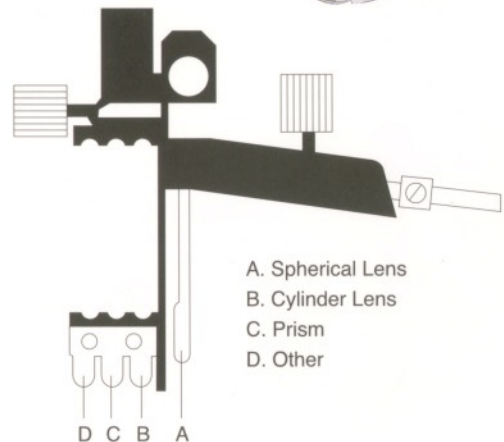

SPECIFICATIONS:

Length of Temple	Min. 100mm / Max. 137mm
PD Adjustment	Min. 100mm / Max. 137mm (Calibrated 24mm to 37mm each side in 1mm steps)
Elevation Distance Adjustment Range	15mm
Vertex Distance Adjustment Range	18mm (Calibrated ± 5 mm by 1mm steps from 12mm vertex)
Weight	76.5g

*Subject to change in design and/or specifications without advanced notice.

IMPORTANT

In order to obtain the best results with this instrument, please be sure to review all user instructions prior to operation.



TOPCON CORPORATION



Certificate No. JP00/18661QA
TOPCON CORPORATION



Certificate No. EC97J1081
TOPCON CORPORATION



TOPCON CORPORATION

75-1 Hasunuma-cho, Itabashi-ku, Tokyo 174-8580, Japan. Phone: 3-3558-2523/2522 Fax: 3-3960-4214 www.topcon.co.jp

TOPCON MEDICAL SYSTEMS, INC.
37 West Century Road, Paramus, NJ 07652, U.S.A.
Phone: 201-261-9450 Fax: 201-387-2710 www.topcon.com

TOPCON CANADA INC.
110 Provencher Avenue, Boisbriand QC J7G 1N1, CANADA
Phone: 450-430-7771 Fax: 450-430-6457 www.topcon.ca

TOPCON EUROPE B.V.
Essebaan 11, 2908 LJ Capelle a/d IJssel, THE NETHERLANDS
Phone: 010-4585077 Fax: 010-4585045 www.topconeuropa.com

ITALY OFFICE
Via Monfalcone 39, 20092 Cinisello B. mo(MI), ITALY
Phone: 02-61-25-583 Fax: 02-61-25-927 www.topcon.it

DENMARK OFFICE
Praestemarkavej 25, 4000 Roskilde, DENMARK
Phone: 046-32-7500 Fax: 046-32-7555

IRELAND OFFICE
Unit 69 Western Parkway Business Centre,
Lower Ballymount Road, Dublin 12, IRELAND
Phone: 01460-0021 Fax: 01460-0129

TOPCON S.A.R.L.
89, Rue de Paris 92585 Clichy, Cedex, FRANCE
Phone: 01-41069494 Fax: 01-47390251

TOPCON DEUTSCHLAND G.m.b.H.
Giesserallee 31-33, D-47877 Willich, Germany
Phone: 02154-8850 Fax: 02154-885111 www.topcon.de

TOPCON SCANDINAVIA A.B.
Neogatan 2 S-43151 Molndal, SWEDEN
Phone: 31-7109200 Fax: 31-7109249

TOPCON ESPAÑA S.A.
Frederic Mompou 5, ED. Euro 3, 08960 Sant Just Desvern Barcelona, SPAIN
Phone: 34-93-4734057 Fax: 34-93-4733932 www.topcon.es

TOPCON (GREAT BRITAIN) LTD.
Topcon House, Kermet Side, Bone Lane, Newbury,
Berkshire RG14 5PX United Kingdom
Phone: 01635-551120 Fax: 01635-551170 www.topcon.com.au

TOPCON SOUTH ASIA PTE.LTD.
Blk 192 Pandan Loop, 407-01 Pantech Industrial Complex, Singapore 128381
Phone: 62780222 Fax: 62733540 www.topcon.com.sg

TOPCON AUSTRALIA PTY.LTD.
Unit 18, 4 Avenue of Americas Newington NSW 2127 AUSTRALIA
Phone: 02-8748-8777 Fax: 02-9647-2926 www.topcon.com.au

TOPCON INSTRUMENTS (MALAYSIA) SDN. BHD.
Exalia Business Park Block C, 1st Floor, Jalan Ampang Putra,
Taman Ampang Hill, 55100 Kuala Lumpur, MALAYSIA
Phone: 03-42701192 Fax: 03-42704508

TOPCON INSTRUMENTS (THAILAND) CO.,LTD.
77/162 Sinn Sathorn Tower, 37th Floor, Krungthorburri Road,
Klongtongnai, Klongsarn, Bangkok 10600 THAILAND
Phone: 02-440-1152-7 Fax: 02-440-1158

TOPCON KOREA CORPORATION
2F Yooseung Bldg., 1595-3, Seocho-Dong, Seocho-Gu, Seoul, 137-876, Korea
Phone: 2-2055-0321-7 Fax: 2-2055-0319 www.topcon.co.kr

TOPCON OPTICAL (H.K.) LTD.
2/F., Meeco Industrial Bldg., No.53-55 Au Pui Wan Street,
Fo Tan Road, Shatin, N.T., Hong Kong
Phone: 852-2690-1328 Fax: 852-2690-2221 www.topcon.com.hk

TOPCON CORPORATION BEIJING OFFICE

1070 Poly Plaza Building, 14 Dongzhimen
Nandajie, Dongcheng District, Beijing, 100027 CHINA
Phone: 10-6501-4191 Fax: 10-6501-4190

TOPCON CORPORATION SHANGHAI OFFICE
Room 901-903 Lippo Plaza No.222 Huai Hai Zhong Lu Shanghai, 20021 CHINA
Phone: 21-62701877 Fax: 21-53966485

TOPCON CORPORATION BEIRUT OFFICE
P.O. BOX 70-1002 Antelias, BEIRUT-LEBANON
Phone: 961-4-523525/523526 Fax: 961-4-521119

TOPCON CORPORATION DUBAI OFFICE
C/O Atlas Medical FZCO, P.O. BOX 54304, C-25,
Dubai Airport Free Zone, United Arab Emirates
Phone: 971-4-2995900 Fax: 971-4-2995901