

QuickSee®

Accurate autorefractor anywhere

More than 2 million people tested in over 30 countries



World's most accurate handheld autorefractor

Binocular and open view

Eliminates patient accommodation for maximum reliability

Wavefront aberrometry

The most comprehensive method to measure ocular refractive errors

Dynamic measurements

Produces results with high confidence

APPLICATIONS

Accurate for the clinic, durable for the field



Mobile eyecare

Home visits • School and office screenings • Nursing homes



Global health

Mission trips • Population health • Initiatives in low resource settings



Eyecare clinics

Patient screening • Post-op • During pretest or refraction

Binocular measurements in 10 seconds

- Accelerates subjective refraction with a more accurate starting point
- Enables high throughput autorefractor

ADA compliant & easy to use

- Ideal for patients with mobility disabilities and/or physical challenges
- Easy to learn
- Works anywhere
- Patient friendly

Field durable

- Calibration free
- Operates in humid and dusty settings
- Can be used indoors and outdoors, in most light settings
- Includes hardened carrying case
- Operates up to 8 hours on battery

QuickSee®

Accurate autorefractometry anywhere

CLINICAL PERFORMANCE SPECIFICATIONS

Intended patient population	5–85 years old
Accuracy (agreement with subjective refraction)	<= 0.25 D: 60-70% of patients <= 0.5 D: 80-90% of patients (see publications)
Cylindrical range	-6D to +6D, increments of 0.01D, 0.125D, 0.25D
Axial range	0–180°, increments of 1, 5, 10 degrees
Spherical range	-10D to +10D, increments of 0.01D, 0.125D, 0.25D
Interpupillary distance range	47–78mm continuous
Accommodation control	Binocular open view
Pupil size	2–8mm
Cycloplegia requirement	None
Amblyopia or strabismus	May require monocular measurement
Contact lens wearers	Over-refraction compatible
Dilation requirement	None
Illumination requirements	Works in any illumination
Furniture requirements	None
Materials	This product is latex-free
Storage/transport	Travel hardshell case included

TECHNICAL SPECIFICATIONS

Wavefront Refraction Engine™ technology patented by PlenOptika	
Technology	Wavefront aberrometry
Measurement mode(s)	Binocular and monocular
Measurement time	10 seconds
Calibration	Factory calibrated, no field calibration needed
Battery life	6–8 hours of operation; 6 hours charge time (10 hours for initial charge)
Communications	Bluetooth via QuickSee Companion App and Bluetooth printer
Measurement capacity	10,000 measurements
Regulatory	Class I product, FDA Class IIa product, CE
Laser safety	Class 1, enclosed
Electrical safety	Compliant with IEC 60601-1(2005, 3.1 ed.)
Weight	3.14lbs / 1.42kg (without case) 6.95lbs / 3.15kg (with case)
Dimensions	6.5 x 11 x 3.25" / 16.5 x 28 x 8.25 cm
Warranty	1 year warranty

Contact


+1 617-862-2203
info@plenoptika.com

United States

PlenOptika, Inc.
955 Massachusetts Avenue
Number 339
Cambridge, MA 02139

Europe

PlenOptika Europe S.L.U.
Parque Científico de Madrid
Calle Faraday 7, Number 007
Madrid 28049, Spain

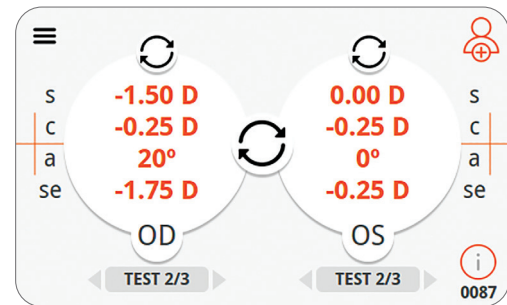

User interface: Simple and intuitive

KEY PUBLICATIONS

1. Durr NJ, Dave SR, Lim D, et al "Quality of eyeglass prescriptions from a low-cost wavefront autorefractor evaluated in rural India: results of a 708-participant field study" British Medical Journal Open Ophthalmology 2019;4:e000225. doi: 10.1136/bmjophth-2018-000225
2. Rubio M, Hernandez CS, Seco E, Perez-Merino P, Casares I, R. Dave SR, Lim D, Durr NJ, Lage E. "Validation of an affordable handheld wavefront autorefractor". Optometry and Vision Science, 2019 (accepted, in press)