

Longlife™ Project-O-Chart® Manual Projector

The standard for manual projectors.



Longlife[™] **Project-O-Chart**[®] Manual Projector

When it was first introduced, the **Longlife POC** set new standards for manual projectors and is currently the most popular manual projector among eye care professionals around the world.

The **Longlife** delivers performance superior to other manual projection systems. Features include:

- Long Halogen bulblife, lasting ten times longer than incandescent bulbs
- Quick-change socket design replaces the old fashioned "push and twist" bulbs that break in your hands
- A bright, white, uniform illumination, ensuring consistent test results
- Pleasing design complements any office decor
- Superior Reichert-quality optics for bright, high-contrast images
- Simple installation fits most existing mounts


Vectographic Slides

Reichert's Vectographic slides create test situations that simulate normal binocular vision through the use of polarized iodine crystal characters in a binocularly viewed structure. This format facilitates fusion of the images cortically. Vectographic slides provide the ability to conduct a wide variety of tests for both children and adults. These include:

- Astigmatism tests
- Acuity tests
- Muscle balance tests (Fixation disparity)
- Acuity balance tests
- Tests to determine malingering
- Stereopsis tests
- Fixation disparity tests with or without paracentral fusion lock



Longlife POC Slide Library

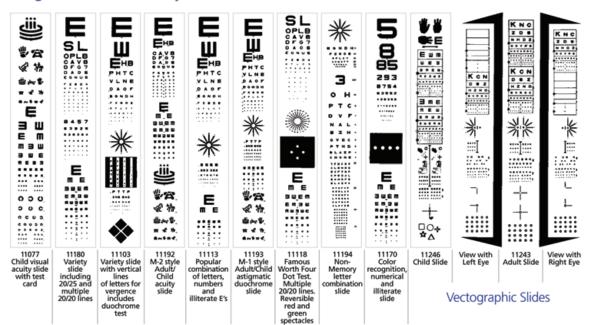

Catalog Numbers

| LONGLIFE POO | : |
|--------------|-----------------------------------|
| 12084 | Longlife POC (head only) |
| 12094 | Yoke (required) |
| 11084 | 13.5" x 13.5" screen |
| 11808 | 20" x 20" non-depolarizing screen |
| 11143 | Halogen Bulb |

Catalog Numbers

| MOUNTS | |
|--------|------------------------|
| 12090 | Table Mount |
| 12091 | Wall Mount |
| 12092 | Instrument Stand Mount |
| 12096 | Floor Stand |



Reichert, Inc.

3362 Walden Avenue, Depew, NY 14043 Toll Free: 888-849-8955 • Tel: (716) 686-4500 Fax: (716) 686-4545 • www.reichert.com

Reichert products are designed and manufactured under quality processes meeting ISO 9001 requirements.